

ZENFLOW™ II HP PTA Balloon Catheter

Catheter Effective Length(cm)			Balloon Diameter [mm]	Balloon Length (mm)	Sheath (F)	NP (ATM)	RBP (ATM)
40	75	135					
ZH030020040	ZH030020075	ZH030020135	3.0	20	5	10	26
ZH040020040	ZH040020075	ZH040020135	4.0		5	10	24
ZH050020040	ZH050020075	ZH050020135	5.0		5	10	24
ZH060020040	ZH060020075	ZH060020135	6.0		5	10	24
ZH070020040	ZH070020075	ZH070020135	7.0		6	10	20
ZH080020040	ZH080020075	ZH080020135	8.0		6	10	20
ZH090020040	ZH090020075	ZH090020135	9.0		7	8	18
ZH100020040	ZH100020075	ZH100020135	10.0		7	8	14
ZH120020040	ZH120020075	ZH120020135	12.0		8	8	14
ZH030030040	ZH030030075	ZH030030135	3.0	30	5	10	26
ZH040030040	ZH040030075	ZH040030135	4.0		5	10	24
ZH050030040	ZH050030075	ZH050030135	5.0		5	10	24
ZH060030040	ZH060030075	ZH060030135	6.0		5	10	24
ZH070030040	ZH070030075	ZH070030135	7.0		6	10	20
ZH080030040	ZH080030075	ZH080030135	8.0		6	10	20
ZH090030040	ZH090030075	ZH090030135	9.0		7	8	18
ZH100030040	ZH100030075	ZH100030135	10.0		7	8	14
ZH120030040	ZH120030075	ZH120030135	12.0		8	8	14
ZH030040040	ZH030040075	ZH030040135	3.0	40	5	10	26
ZH040040040	ZH040040075	ZH040040135	4.0		5	10	24
ZH050040040	ZH050040075	ZH050040135	5.0		5	10	24
ZH060040040	ZH060040075	ZH060040135	6.0		5	10	24
ZH070040040	ZH070040075	ZH070040135	7.0		6	10	20
ZH080040040	ZH080040075	ZH080040135	8.0		6	10	20
ZH090040040	ZH090040075	ZH090040135	9.0		7	8	18
ZH100040040	ZH100040075	ZH100040135	10.0		7	8	14
ZH120040040	ZH120040075	ZH120040135	12.0		8	8	14
ZH030060040	ZH030060075	ZH030060135	3.0	60	5	10	26
ZH040060040	ZH040060075	ZH040060135	4.0		5	10	24
ZH050060040	ZH050060075	ZH050060135	5.0		5	10	24
ZH060060040	ZH060060075	ZH060060135	6.0		5	10	24
ZH070060040	ZH070060075	ZH070060135	7.0		6	10	20
ZH080060040	ZH080060075	ZH080060135	8.0		6	10	20
ZH090060040	ZH090060075	ZH090060135	9.0		7	8	18
ZH100060040	ZH100060075	ZH100060135	10.0		7	8	14
ZH120060040	ZH120060075	ZH120060135	12.0		8	8	14
ZH030080040	ZH030080075	ZH030080135	3.0	80	5	10	26
ZH040080040	ZH040080075	ZH040080135	4.0		5	10	24
ZH050080040	ZH050080075	ZH050080135	5.0		5	10	24
ZH060080040	ZH060080075	ZH060080135	6.0		5	10	24
ZH070080040	ZH070080075	ZH070080135	7.0		6	10	20
ZH080080040	ZH080080075	ZH080080135	8.0		6	10	20
ZH090080040	ZH090080075	ZH090080135	9.0		7	8	18
ZH100080040	ZH100080075	ZH100080135	10.0		7	8	14
ZH120080040	ZH120080075	ZH120080135	12.0		8	8	14
ZH030100040	ZH030100075	ZH030100135	3.0	100	5	10	26
ZH040100040	ZH040100075	ZH040100135	4.0		5	10	24
ZH050100040	ZH050100075	ZH050100135	5.0		5	10	24
ZH060100040	ZH060100075	ZH060100135	6.0		6	10	24
ZH070100040	ZH070100075	ZH070100135	7.0		6	10	20
ZH080100040	ZH080100075	ZH080100135	8.0		6	10	20
ZH030120040	ZH030120075	ZH030120135	3.0	120	5	10	26
ZH040120040	ZH040120075	ZH040120135	4.0		5	10	24
ZH050120040	ZH050120075	ZH050120135	5.0		5	10	24
ZH060120040	ZH060120075	ZH060120135	6.0		6	10	22
ZH070120040	ZH070120075	ZH070120135	7.0		6	10	18
ZH080120040	ZH080120075	ZH080120135	8.0		6	10	18
ZH030150040	ZH030150075	ZH030150135	3.0	150	5	10	26
ZH040150040	ZH040150075	ZH040150135	4.0		5	10	24
ZH050150040	ZH050150075	ZH050150135	5.0		5	10	24
ZH060150040	ZH060150075	ZH060150135	6.0		6	10	22
ZH070150040	ZH070150075	ZH070150135	7.0		6	10	18
ZH080150040	ZH080150075	ZH080150135	8.0		6	10	18
ZH030180040	ZH030180075	ZH030180135	3.0	180	5	10	26
ZH040180040	ZH040180075	ZH040180135	4.0		5	10	24
ZH050180040	ZH050180075	ZH050180135	5.0		5	10	24
ZH060180040	ZH060180075	ZH060180135	6.0		6	10	22
ZH070180040	ZH070180075	ZH070180135	7.0		6	10	18
ZH080180040	ZH080180075	ZH080180135	8.0		6	10	18
ZH030200040	ZH030200075	ZH030200135	3.0	200	5	10	26
ZH040200040	ZH040200075	ZH040200135	4.0		5	10	24
ZH050200040	ZH050200075	ZH050200135	5.0		5	10	24
ZH060200040	ZH060200075	ZH060200135	6.0		6	10	22
ZH070200040	ZH070200075	ZH070200135	7.0		6	10	18
ZH080200040	ZH080200075	ZH080200135	8.0		6	10	18

ZENFLOW™ II HP PTA Balloon Catheter is not available in markets requiring FDA clearance and CE certification.

Zylox-Tonbridge Medical Technology Co., Ltd

Zylox-Tonbridge Industrial Park, No. 270 Shuyun Road, Yuhang District, Hangzhou, Zhejiang Province, China.

TEL: (+86) 571 8861 0082

FAX: (+86) 571 8861 0608


E-mail: info@zyloxmedical.com

<https://www.zyloxtb.com/>

©2024 Zylox-Tonbridge Medical Technology Co.,Ltd



ZENFLOW™ II HP
PTA Balloon Catheter

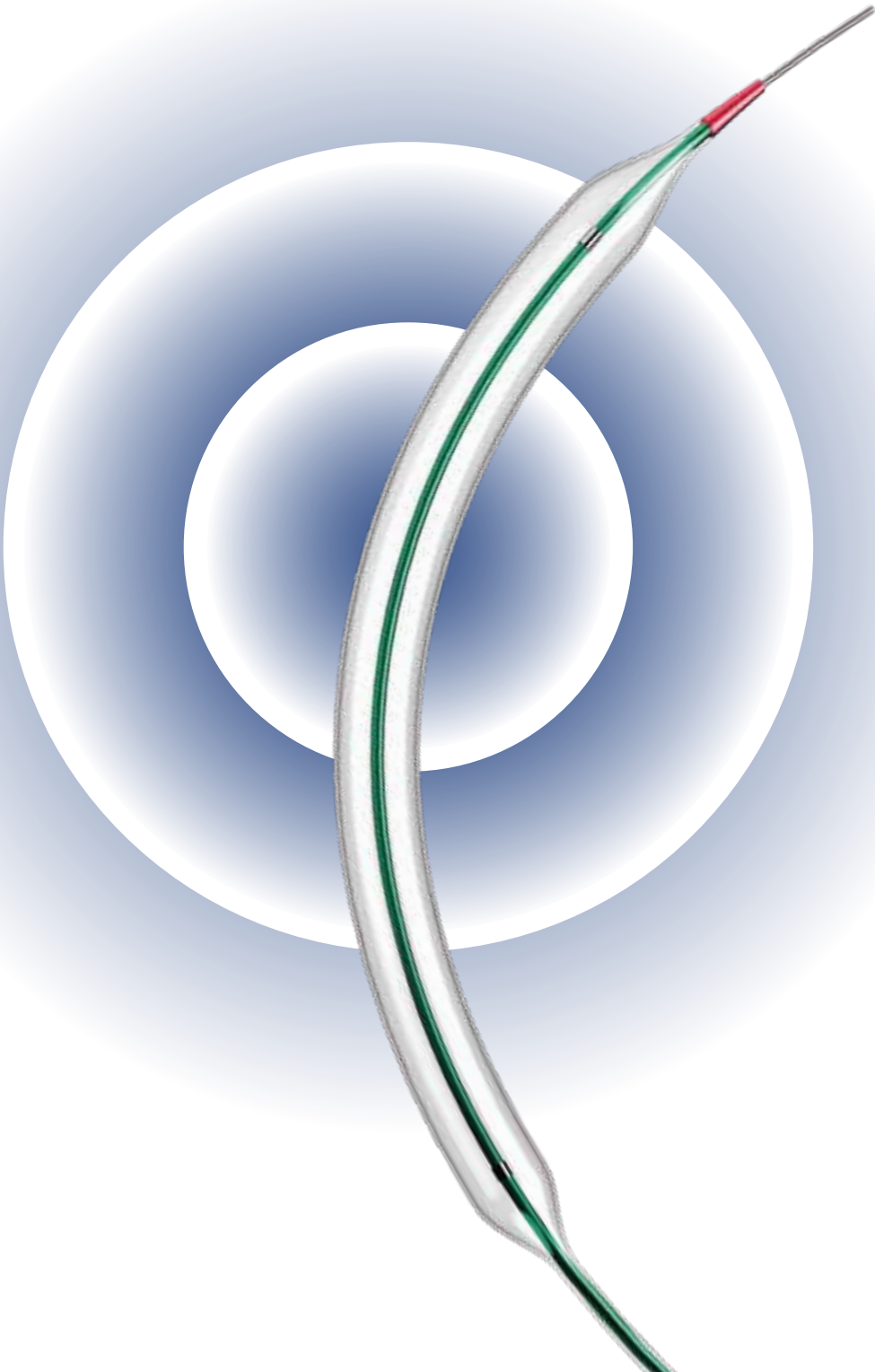

ZENFLOW™ II HP

PTA Balloon Catheter

Powerful Solution For Resistant Lesions!

Up to **26**
ATM RBP



6-Folding



Dual Lumen
Design



Deliver stronger outcomes for both AV and vascular access patency

- **Balloon Materials** (Nylon + Pebax)
Optimized Wall Thickness: The Nylon and Pebax composition allows for a thinner wall while maintaining high RBP and low compliance.
- **Advanced Folding Technology**
Consistent Low Profile: Maintains shape for easy insertion and navigation, ensuring efficient rewinding after deflation.
- **Soft tip; low entry profile**
Easy and Safe Access: Allows atraumatic entry into complex structures, minimizing the risk of vessel injury.
- **Dual Lumen Shaft Design**
Enhanced Performance: Enables fast deflation and offers improved push-ability during advancement.
- **Wide range of sizes**
 - With a maximum diameter of 12 mm, and up to 200 mm balloon length.
 - Available in 40 cm, 75 cm, and 135 cm shaft lengths

Technology data

Catheter Type: Over-The-Wire (OTW)

Guidewire: 0.035"

Folding: Three/six wings

Balloon Coating: Hydrophilic coating

Balloon Markers: 2